



CareLink™ Follow-up

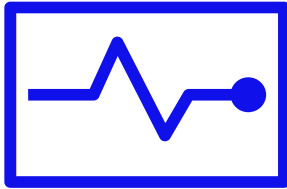

9

Goals

HCP quick reference

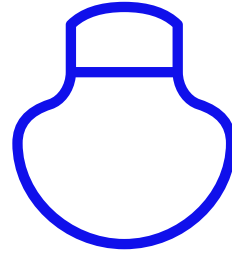

1

Report

Assessment & Progress



1

Question

Are they meeting goals?



Metric	Time above ² 13.9 mmol/L	Time above ² 10.0 mmol/L	Time in Range ² 3.9 - 10.0 mmol/L	Time below ² 3.9 mmol/L	Time below ² 3.0 mmol/L	A1C ³	Coefficient of variation (CV) ⁴	SmartGuard™ Use	Sensor Use ⁵
Goals	< 5%	< 25%	70%	< 4%	< 1%	< 7.0% Adults (< 53 mmol/mol)	< 36%	≥ 85%	≥ 85%

CareLink™ Follow-up



1

Question

Are they meeting goals?

YES?

You're done - yay!

NO?

1. Is the system working as hard as it can be?

SmartGuard™ target @ 5.5mmol/L

Active Insulin Time @ 2h

Bolus Increment @ 0.025

2. Are carb ratios optimized?